

™ Simple
Connected
Powerful




MTS101

Monitor with Touch Screen

The MTS101 is a medical grade 10.1" Multi-Touch Screen monitor, which is USB-powered and can be used in the operating rooms (OR) for controlling any MediCapture MVR or other OEM imaging systems.

Designed for modern operating rooms and surgery centers, the vendor-neutral MTS101 enhances the user experience by providing comfortable monitor positioning and a larger user interface for entering patient data, editing videos and developing reports.

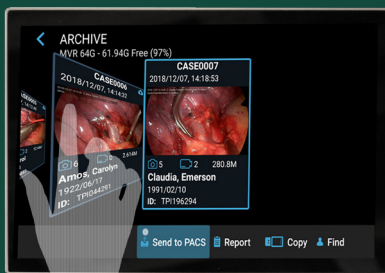

■ Simple to Use

With MTS101 and its 5-points multi-touch, you can handle your workflow with a swipe of finger -- just like your smartphone.

■ Designed for the OR


The MTS101 eliminates operational boundaries inside the operating rooms. With its comfortable 10.1" (26cm) multi-touch screen and easy VESA Mount capability, you can choose any convenient location inside the OR and control any imaging system.

■ Easy to connect & Eco-friendly

Powered by USB Type-C, the MTS101 consumes less energy than traditional monitors, supporting you in reducing overall energy consumption and lowering operating costs -- Just Plug & Play.

Medical Imaging Made Easy™

MTS101 Features

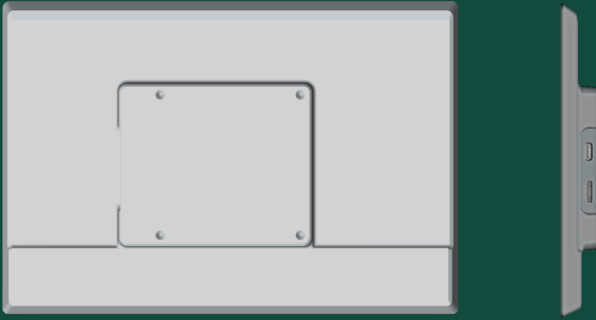

Backside with VESA 75

Side connections



Medical grade

It provides accurate quality, safety, intuitive user interaction, and easy-to-clean surfaces that help maintain a safe and hygienic environment in OR.



Vendor-neutral

Using MTS101 with all the OEM-specific equipments in OR can help streamline workflows and minimize compatibility issues.



Flexible to mount in OR

With VESA 75, various mounting options are available, such as support arms, table stands, and clamps to poles or rails³.

MTS101 Specifications¹

Device Type	Medical Monitor with Touch Screen Aluminium Alloy, 26cm / 10.1" diagonal, 260 cd/m ² , Aspect Ratio: 16:10, Multi-Touch 5 points
Monitor Resolution	1920 x 1200
Video Inputs	Mini HDMI (Type-C) Supports HDMI (DVI via adapter)
Input Resolution	1920x1200p60
USB	Type-C 2.0 for power and control Compatible to USB 3.0 host
Mechanical	Size: 244,1 x 167,1 x 20,2 mm / 9,6 x 6,6 x 0,8" Weight: 0,63 kg (1,4lbs)
Environment	Storage: -40° to +85° C Operating: -20°C to +40°C
Power	From USB Type-C, DC 5V / 0.8A
Certifications	EN 60601-1:2006 + A1:2013 + A12:2014 - Medical Equipment Safety, UL, MDR 2017/745/EU, Class 1, EN 60601-1-2:2015 - Electromagnetic Compatibility EMC, FCC Part 15B - Class B, RoHS, WEEE, UKCA
Order Number	MTS101 : MTS monitor, HDMI to Mini HDMI cable, USB Type-A to USB Type-C cable, Information for Use (English) ² ST004MTS : Acryl Stand For MTS101



¹ Subject to change without notice

² For more languages of the User Guide visit www.medicapture.com/support

³ MediCapture exclusively sells the Acryl Stand for mounting purposes and does not sell any other VESA 75 mounting options.

MediCapture, Inc. – North America

2250 Hickory Road, Suite 200
Plymouth Meeting, PA 19462, USA

+1 888 922 7887

www.medicapture.com

sales@medicapture.com

MediCapture Europe B.V.

Stationsplein 5
6131 AT Sittard, The Netherlands

+31 46 208 1120

www.medicapture.com

sales.europe@medicapture.com

MediCapture APAC Inc.

3FL-6, No. 432, Section 1, Keelung Road,
Xinyi District 110, Taipei City, Taiwan, Republic Of China

+886 2 2345 3668

www.medicapture.com

sales.asiapacific@medicapture.com